

The Standard for Pathology Imaging

Helping pathology partners develop roadmaps

- Professional images from streamlined workflows
- Shared central image archive
- Secure telepathology interface
- Synchronized LIS connectivity

Based on the principles of:

- Practical, budget positive solutions that pay for themselves
- Integrated, stepwise implementation
- Fitting your existing infrastructure of equipment and software
- Solid performance 55+ years of problem-solving experience



SPOT Imaging Partners



QR to SPOT products



6540 Burroughs Avenue • Sterling Heights • MI • 48314-2133
586.731.6000 • toll free (US) 1-866-604-SPOT • spotimaging.com • info@spotimaging.com

PT 8PG BRO. 05.24



**SAVES TIME
& ELIMINATES
MISTAKES!**

PathTracker™ Bulk Barcode Reader



IN PROCESS SAMPLE TRACKING
VISUAL CONFIRMATION OF BARCODE SCANS
LOGS BARCODES INTO YOUR LIS

Targeted Solutions

PathTracker™

Hand Scanning Batches of Barcodes is Problematic...

Every lab manager and pathology director knows the importance of tracking samples from start to finish. A few losses and your lab's credibility will be in question from both your clinical colleagues as well as your state's regulatory agency.

Best practices dictate "way point" scanning of samples as they traverse histological processing, pathology diagnosis and long term archiving.

While "way point" scanning has helped sample tracking greatly, it has introduced its own issues

related to the number of scans required at each process step and the difficulty of ensuring all cassettes and slides in a basket or index are scanned. The result is mind numbing work and frustrating re-scanning of mismatched manifests.



PathTracker automates this process, eliminating non-value added work and its associated errors.

Bulk Sample Tracking Made Easy


PathTracker replaces laborious and error prone hand scanning of cassette and slides with a single push of a button. PathTracker incorporates patented technology to acquire, process, analyze and log all the barcodes in the field of view.

150 cassette processing baskets are scanned in less than 9 seconds, compared to 7.5 minutes for hand scanning. Slide indexes are completed in less than 2 seconds vs 1 minute 20 seconds for hand scanning.


Critically, you can rest assured that all the cassettes were read and any damaged or poorly printed barcodes are flagged. PathTracker then provides a set of correction tools to automatically or manually correct these barcodes ensuring your production run continues through its processing.

Once approved, PathTracker transfers the data set to your LIS and stores the source image and data set to its backup archive.

Save Time

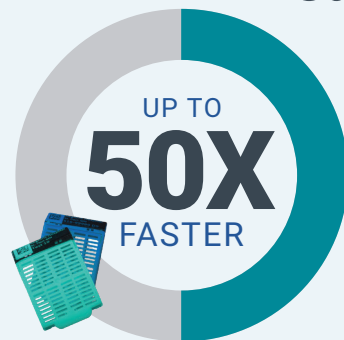

150 CASSETTE BASKET
PATHTRACKER = < 9 Seconds
HAND SCANNING = > 450 Seconds

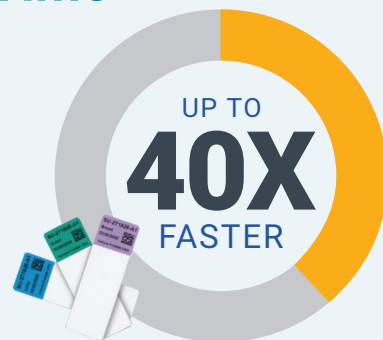

20 SLIDE FOLDER
PATHTRACKER = < 2 Seconds
HAND SCANNING = > 80 Seconds

Embrace the Digital Pathology Promise




- Minimize Errors
- Consistently Track Samples: Origin to Archive
- Increase Productivity and Lab Utilization
- Eliminate Operator Frustration
- Maintain Personnel Focus

PathTracker™ | Hardware

Specifications

SPECIFICATIONS	METRIC	SAE
Computer OS	Win 10 Pro	Win 10 Pro
USB 3.0 Ports	4	4
LAN Port	1 GigaBit Cat5	1 GigaBit Cat5
Washable Keyboard and Mouse	(Standard)	(Standard)
Hand Held Barcode Reader	(Standard)	(Standard)
DISPLAY		
Touch Screen	(Standard)	(Standard)
Size (diagonal)	54.5cm	21.5"
Resolution	1920 x 1080	1920 x 1080
SAMPLE CAPACITY		
Maximum Cassette FOV (HxV)	22.9 x 15.3 cm	9.0" x 6.0"
Maximum Number of Cassettes	150	150
Maximum Microscope FOV (HxV)	34.3 x 22.9 cm	13.5" x 9.0"
Maximum Number of Slides	20	20
Maximum Container Height	14.0 cm	5.5"
INSTRUMENT SIZE		
Height	81.3 cm	32.0"
Width	52.3 cm	20.6"
Depth	61.0 cm	24.0"
WEIGHT		
Instrument	38.2 Kgs	84 lbs
Shipping	88.4 Kgs	195 lbs
VENTILATION COUPLING		
Diameter	10.2 cm	4"
Height from Table Top	21.3 CM	8.4"
Weight	38.21 Kgs	84.25 lbs
POWER		
Input	AC100~240V 6A 50/60Hz	
Fuse	8 AMP / 240V Slow Blow, 5x20mm	
Operation Temperature	0-40°C or 32°- 104° F	
Caution	For indoor use only. For use with Diagnostic Instruments Inc. PathTracker Software only.	
Certifications	CE, FCC Class A, RoHS	

Dimensions

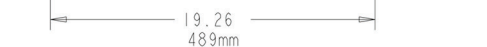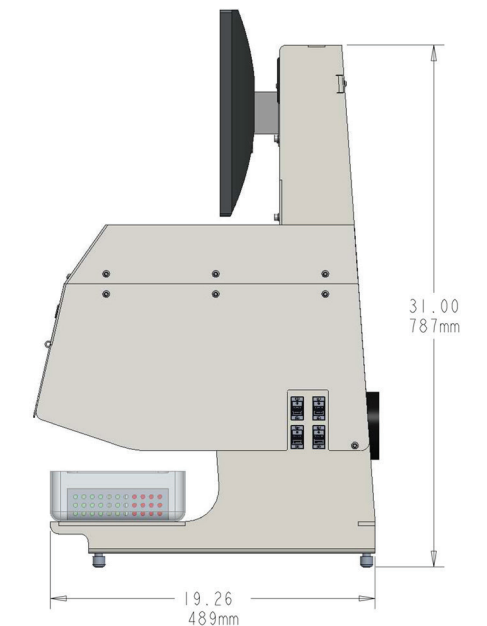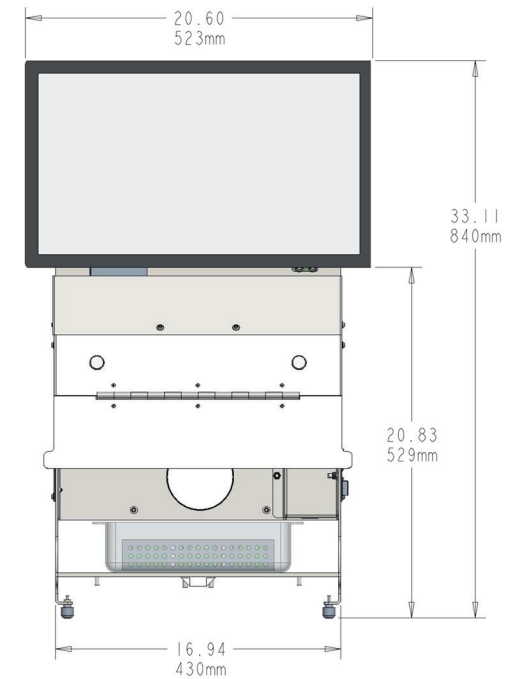

PathTracker™ | Key Value Generating Features

Return on Your Investment... Fast Scans Save Time and Money!

Save Time Save Money

For 120,000 Cassettes

SAVE:
567 EFT Hours¹
Save \$34,000³

For 600,000 Slides

SAVE:
333 EFT Hours²
Save \$20,000³

PathTracker Return on investment starts with providing fast scans that allow operators to complete all their check out or check in scans in the time they would have spent checking only 1 bulk container. This savings can free up your personnel to perform more meaningful work, and eliminate overtime.

Secondly, **PathTracker** is so easy to use that it eliminates the common lab practice of having "The Person" who is qualified careful enough to perform the bulk scanning for the lab. **PathTracker** is one more step on your journey to lab efficiency and sample accountability.

Reduce Errors & Speed Work

PathTracker Leverages Simple Facts:

FACT 1: There is a baseline human error rate. The reduction of manual operations reduces errors...
PathTracker bulk scanning reduces error frequency by:

- A factor of 150:1 for cassettes
- A factor of 20:1 for slides


FACT 2: Human operators can identify errors more easily and quickly in a graphical presentation verses a textual list of codes...

- PathTracker presents its scan results in a visual presentation with markers to indicate recorded barcodes and barcodes to be corrected.
- If a scanning error does happen, it can be quickly identified and addressed without having to rescan all the samples in a container, saving time and keeping your team on schedule.

FACT 3: The sooner a sample is flagged as late or missing, the more limited the search area will be...

- PathTracker provides notification of missing or late samples quickly, in process, bringing to attention errant sample container as they occur rather than days later at the end of the process.
- When a sample is lost It is easy to find with PathTracker's scan attributes that quickly identify operator, courier name, processor used and processor program.

Reduce Human Error

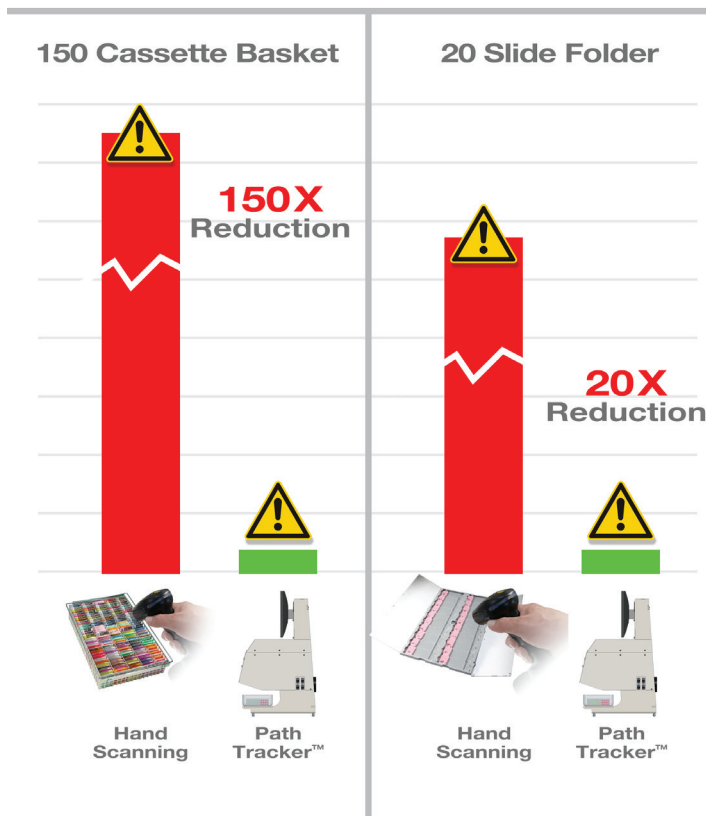

¹ Carraro P, Plebani M. Errors in stat laboratory: types and frequencies 10 years later. Clin Chem 2007; 53: 1338-42.
² Effectiveness of Barcode Tracking in Documenting Errors & Preventing Patient Specimen Identification Incidents in Anatomical Pathology Laboratory, Dr .K.C. Lee, Department of Pathology Princess Margaret Hospital, 15 May 2013 PPT
³ Based on 3 seconds per scan, 6 scanning way points

PathTracker™

Archiving Made Manageable | PathTracker


PathTracker is a practical and affordable management system for your slide and cassette archive. First, it saves time scanning your cassettes and slides into the archive system. Next it saves money with the ability to use standard cassettes storage drawers rather than expensive proprietary storage drawers that will be archived for years to come. When you need a cassette or slide, **PathTracker** provides an easy to use search interface that provides the location of the sample for retrieval and then tracks the location of the loan and when its due back. When it's time to return the block or slide, **PathTracker** eliminates the refiling fees. Simply add the slide or cassette to the next storage container and **PathTracker** logs it in, updates the database and tracks its new location.



Backed Up and Available Network Wide

PathTracker is configurable to be a stand-alone instrument or a networked resource that utilizes the IT infrastructure. As a network device, **PathTracker** supports a centralized database for tracking sample assets from multiple terminals across your network. This ensures data backup security as well as remote access from any **PathSuite** system in the network. Automated archive maintenance is configurable with settings for storage duration on screen images, scannable images, tracking data and long term sample transaction data. **PathTracker** is an integrated pathology data solution.

Implementation Is Easy

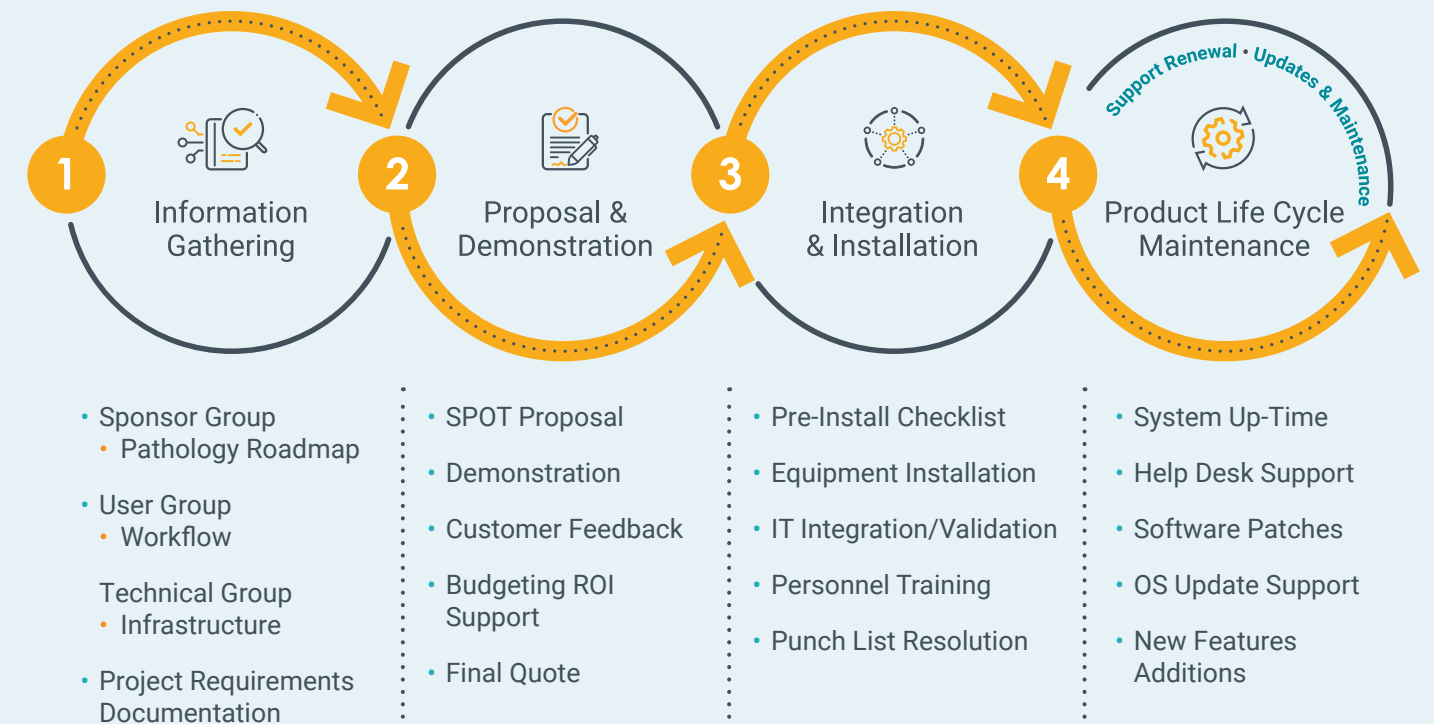

PathTracker™ | Software


Make Bulk Scanning Easy

PathTracker performs its scanning by imaging your bulk sample container, scanning the field for barcodes, and recording the barcodes to the PathTracker and LIS database archives.

Critically, PathTracker's visual interfaces shows you an image of the scanned samples with all the scanned barcodes marked. Unscanned cassettes are easily identified and quickly addressed using PathTracker's correction tools ensuring your production run continues with its processing. During reconciliation scans, immediate notification is provided for missing samples ensuring issues are addressed in timely fashion.

Once a scan is approved, PathTracker transfers the scan data to the PathTracker database as well as to your LIS barcode tracking system.

PathTracker is Fast, Easy and Accurate!

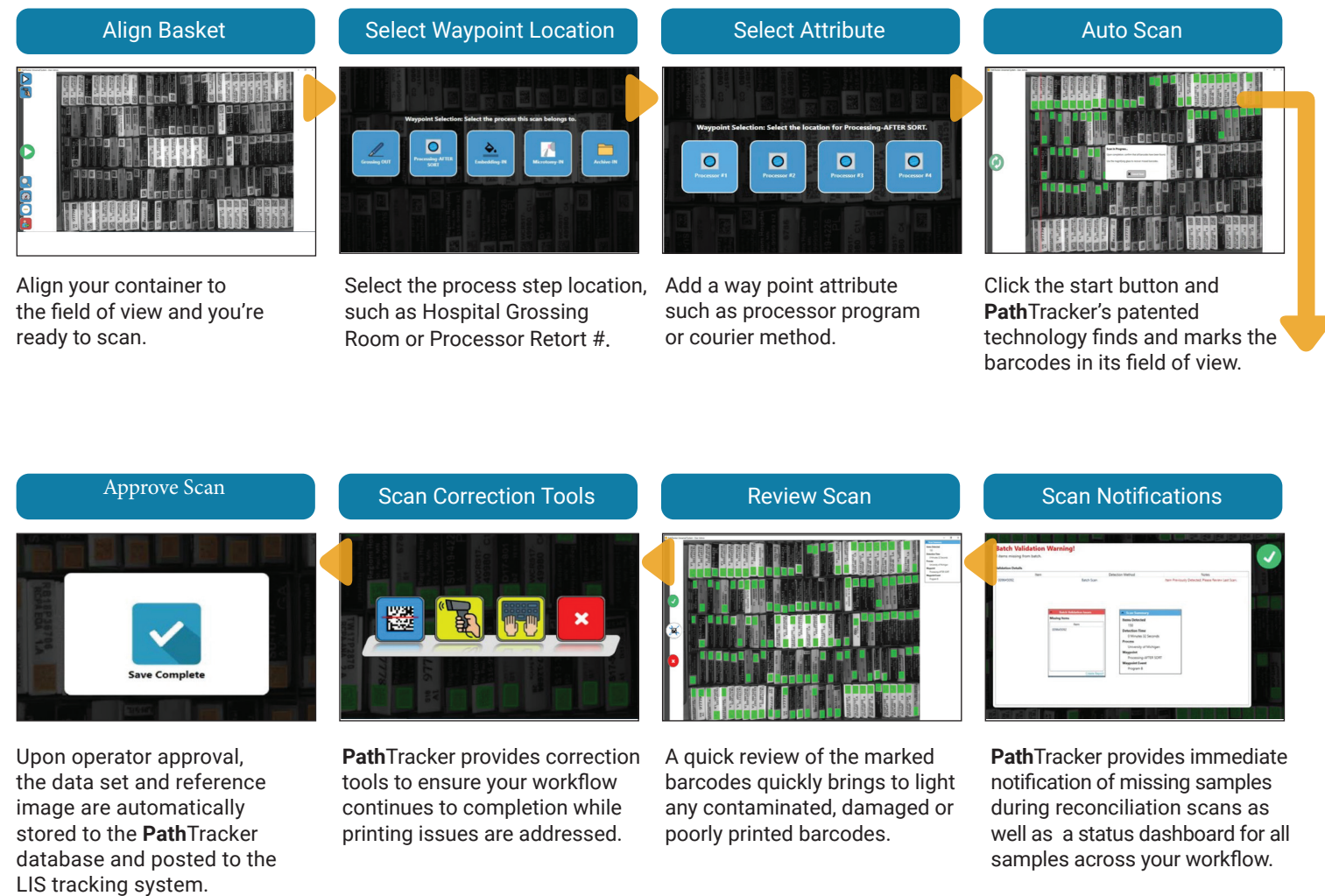


PathTracker™ | Posts Scan Results to your LIS

PathTracker provides automatic posting of sample barcode scans to your LIS via HL7 Message or REST-API / Web Service interfaces. This ensures that PathTracker is keeping your LIS Barcode tracking system updated with the latest locations of each and every sample. Contact your SPOT Imaging representative for the availability and requirements for your LIS tracking software.


Late & Missing Sample Warnings

PathTracker provides missing sample alerts for direct shipping transfers. Checkpoints that have resorting or splitting of the container batch are monitored with transit timers providing late or missing alerts.

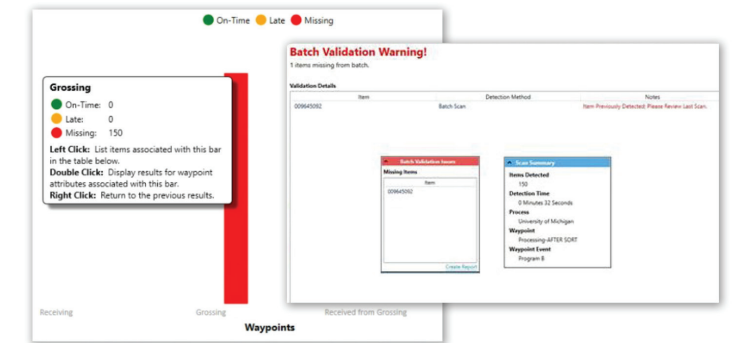

Workflow Modeler

The PathTracker Workflow Modeler maps your institution's Histology processing steps from sample container printing at outlying hospitals through "in lab" tasks and ending with the sample filing room. Once mapped, any step can be easily activated as a tracking way-point with the click of a check box and items added to record attributes such as operator, courier used, or processor program run expanding the data recorded as needed.

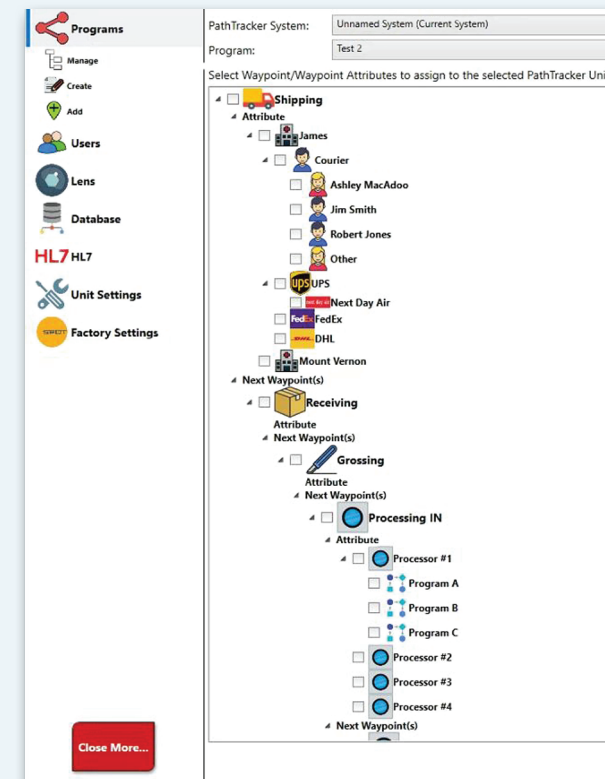

Sample Tracking Management Tools

PathTracker includes a Workflow Dashboard and Report Generator. These tools allow operators to monitor the overall workflow as well as track individual samples as necessary. The report generator tools enable:

- A visual overview of the sample load across the workflow
- Searches for the current status as well as tracking history of individual samples
- Data drill down on a particular waypoint or cohort of samples for workload analysis
- Reporting based on search criteria such as scan date range or waypoints traversed.

